

## FRIENDS Serial Communication GUI

FRIENDS Serial Communication

Status Bar  
Connected

COM  
COM9 - USB Serial Port (COM9) Reload

Baudrate (bps)  
115200 Connect Disconnect

Minimum puffing duration (s) and plot type  
0.3 Stem send

Input  
Read Time Set time Erase Flash Start Data Collection Read Data Yes

RX  
Internal timestamp: 6480AB20 648A  
Number of records: 35  
Timestamps:  
30006478DA670457  
40006478DA674C81  
30006478DA680F26  
40006478DA6866C9  
30006478DA68EC39  
40006478DA692AA2  
10006478DA6C6F6E  
20006478DA6CCEB6  
10006478DA6EAD9  
20006478DA6F024D  
10006478DA7147AA  
20006478DA718AD8  
30006478DA778875  
40006478DA77E029  
30006478DA790309  
40006478DA794AC0  
30006478DA79EA16  
Clear

**Status Bar:** It shows the connectivity status between the computer and the FRIENDS device. It turns into “green” and shows “connected” when the connection between the two devices have been established after selecting the appropriate COM port and baud rate.

**COM:** Appropriate COM port should have been selected for establishing the connection for a particular FRIENDS device.

**Baudrate:** For the FRIENDS device, the baudrate will be 115200.

**Reload:** If the appropriate comp port for the device doesn’t show in the options, reload button should be pressed to get the appropriate comp port.


**Connect:** After setting the comp port and baud rate, connect button should be pressed to establish the connection between the computer and the FRIENDS device

**Disconnect:** To turn off the connection, disconnect button should be pressed.

**Minimum puffing duration:** By default, the value of this entry is 0.0. By setting minimum puffing duration, user can separate the puffing events in the plots based on this threshold.

Example:

For minimum puffing duration 0.1s:

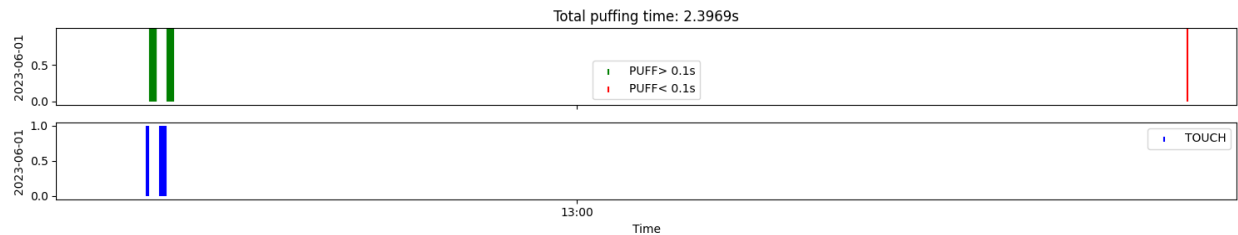

**Plot types:** Three plot types (Stem, Step, and Line) are available for user selection for generating the plots.

**Read Time:** This button returns the device's actual time in POSIX format.

**Set Time:** This button set the local time in POSIX format to the device.

**Erase flash:** This button will erase the device's flash memory. It requires a confirmation (YES) and "Yes" button will be enabled after pressing this button. At the beginning, the erase button will be disabled. After saving the data and generating the plots by "Read Data" button, the erase button will be enabled.

**Start Data Collection:** This button should be pressed before start collecting the data. It will set the local time in POSIX format to the device and erase the flash memory of the device and make it ready to start the data collection.

**Read Data:** This button has multiple tasks. It reads the flash memory, generate 3 text files, and also generate the graphs from the data.


- File 1: It contains the local time in human readable format (e.g., Local Time: 2023-06-05 13:06:47.836938), device's internal timestamp (e.g., Internal timestamp: 647E2437 E3EB), and event's timestamps in **POSIX** format in flash memory.
- File 2: it contains the local time in human readable format (e.g., Local Time: 2023-06-05 13:06:47.836938), device's internal timestamp (e.g., Internal timestamp: 647E2437 E3EB), and event's timestamps in **human readable format** in flash memory.
- File 3: It contains a data table with four columns (Event, Date, Range, Duration\_in\_seconds). One can find each complete event's occurring date, time range and duration from this table. Example:

Event	Date	Range	Duration_in_seconds
Touch	2023-06-01	12:50:31.016953-12:50:31.298843	0.28188999999838416
Touch	2023-06-01	12:50:32.059174-12:50:32.401505	0.34233099999983096
Touch	2023-06-01	12:50:32.922745-12:50:33.166534	0.24378900000010617
PUFF	2023-06-01	12:50:36.435272-12:50:36.807465	0.3721929999956046
PUFF	2023-06-01	12:50:38.667374-12:50:39.008987	0.34161300000414485
PUFF	2023-06-01	12:50:41.279938-12:50:41.542358	0.2624199999991106

After pressing the “Read Data” button, a file dialog will be popped up and it will require a path and a filename to save that 3 files. File 1 will be saved by the given name, file 2 will be saved as given name\_converted and file 3 will be saved as given name\_duration.

After saving the 3 files, the system will generate plots from the data. The number of plots depends on the number of days in the data (one plot for one day). Time from 00:00 to 24:00 represents the horizontal axis x.


Example:

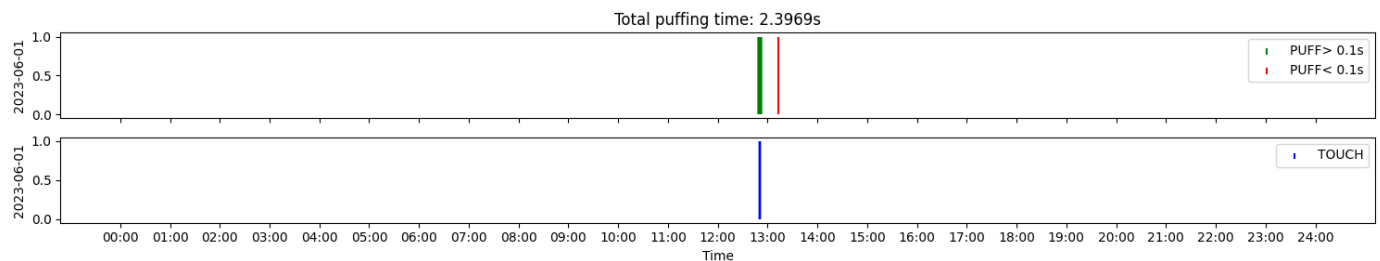

**Clear:** This button will clear all the text in RX monitor

**Save the plot:** For saving the plot, press the icon that is marked by red circle below.

Figure 1

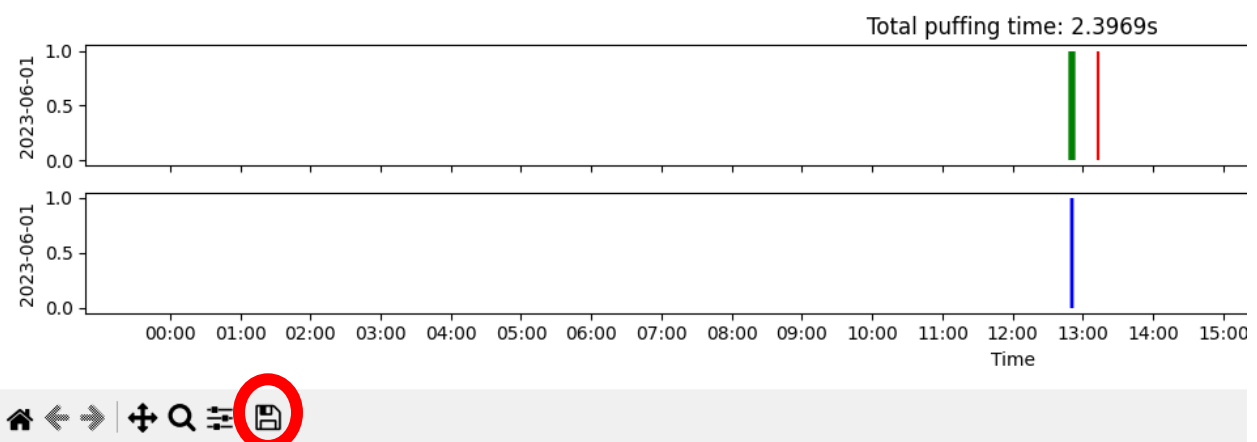

N.B: If any issue occurs while reading the data, please clear the Rx monitor by pressing clear button and try again. This may solve the issue.