

## FRIENDS Serial Communication GUI

The screenshot shows the FRIENDS Serial Communication GUI. The window title is "FRIENDS Serial Communication". At the top, there is a "Status Bar" with a green background and the text "Connected". Below the status bar, there is a "COM" dropdown menu showing "COM9 - USB Serial Port (COM9)" and a "Reload" button. Underneath, there is a "Baudrate (bps)" field set to "115200" with "Connect" and "Disconnect" buttons. Below that, there is a "Minimum puffing duration (s) and plot type" section with a field set to "0.3" and a dropdown menu set to "Stem", with a "send" button. At the bottom, there is an "Input" section with buttons for "Read Time", "Set time", "Erase Flash", "Start Data Collection", "Read Data", and "Yes". The bottom-most section is labeled "RX" and contains a list of internal timestamps: "Internal timestamp: 6480AB20 648A", "Number of records: 35", "Timestamps:", and a list of 16 timestamps. A "Clear" button is located at the bottom right of the RX section.

**Status Bar:** It shows the connectivity status between the computer and the FRIENDS device. It turns into “green” and shows “connected” when the connection between the two devices have been established after selecting the appropriate COM port and baud rate.

**COM:** Appropriate COM port should have been selected for establishing the connection for a particular FRIENDS device. (User needs to install the VCP driver (<https://fdichip.com/drivers/vcp-drivers/>) to work with the FRIENDS device. When FRIENDS is connected to the computer by USB, a serial port appears after installing the VCP driver.)

**Baudrate:** For the FRIENDS device, the baudrate will be 115200.

**Reload:** If the appropriate comp port for the device doesn’t show in the options, reload button should be pressed to get the appropriate comp port.


**Connect:** After setting the comp port and baud rate, connect button should be pressed to establish the connection between the computer and the FRIENDS device

**Disconnect:** To turn off the connection, disconnect button should be pressed.

**Minimum puffing duration:** By default, the value of this entry is 0.0. By setting minimum puffing duration, user can separate the puffing events in the plots based on this threshold.

Example:

For minimum puffing duration 0.1s:

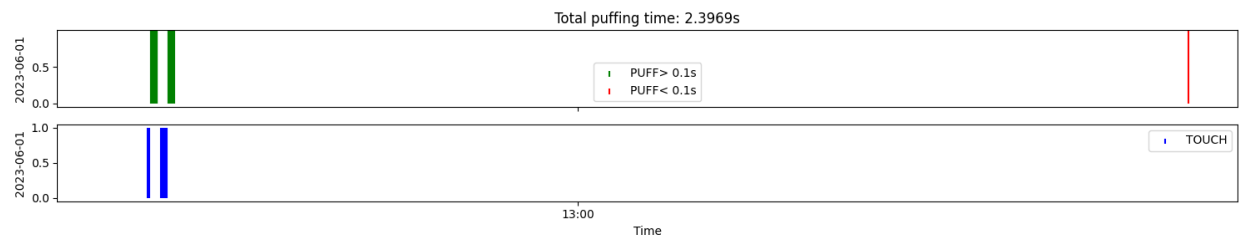

**Plot types:** Three plot types (Stem, Step, and Line) are available for user selection for generating the plots.

**Read Time:** This button returns the device's actual time in POSIX format.

**Set Time:** This button set the local time in POSIX format to the device.

**Erase flash:** This button will erase the device's flash memory. It requires a confirmation (YES) and "Yes" button will be enabled after pressing this button. At the beginning, the erase button will be disabled. After saving the data and generating the plots by "Read Data" button, the erase button will be enabled.

**Start Data Collection:** This button should be pressed before start collecting the data. It will set the local time in POSIX format to the device and erase the flash memory of the device and make it ready to start the data collection.

**Read Data:** This button has multiple tasks. It reads the flash memory, generate 3 text files, and also generate the graphs from the data.


- File 1: It contains the local time in human readable format (e.g., Local Time: 2023-06-05 13:06:47.836938), device's internal timestamp (e.g., Internal timestamp: 647E2437 E3EB), and event's timestamps in **POSIX** format in flash memory.
- File 2: it contains the local time in human readable format (e.g., Local Time: 2023-06-05 13:06:47.836938), device's internal timestamp (e.g., Internal timestamp: 647E2437 E3EB), and event's timestamps in **human readable format** in flash memory.
- File 3: It contains a data table with four columns (Event, Date, Range, Duration\_in\_seconds). One can find each complete event's occurring date, time range and duration from this table. Example:

| Event | Date       | Range                           | Duration_in_seconds |
|-------|------------|---------------------------------|---------------------|
| Touch | 2023-06-01 | 12:50:31.016953-12:50:31.298843 | 0.28188999999838416 |
| Touch | 2023-06-01 | 12:50:32.059174-12:50:32.401505 | 0.34233099999983096 |
| Touch | 2023-06-01 | 12:50:32.922745-12:50:33.166534 | 0.24378900000010617 |
| PUFF  | 2023-06-01 | 12:50:36.435272-12:50:36.807465 | 0.37219299999956046 |
| PUFF  | 2023-06-01 | 12:50:38.667374-12:50:39.008987 | 0.34161300000414485 |
| PUFF  | 2023-06-01 | 12:50:41.279938-12:50:41.542358 | 0.2624199999991106  |

After pressing the “Read Data” button, a file dialog will be popped up and it will require a path and a filename to save that 3 files. File 1 will be saved by the given name, file 2 will be saved as given name\_converted and file 3 will be saved as given name\_duration.

After saving the 3 files, the system will generate plots from the data. The number of plots depends on the number of days in the data (one plot for one day). Time from 00:00 to 24:00 represents the horizontal axis x.


Example:

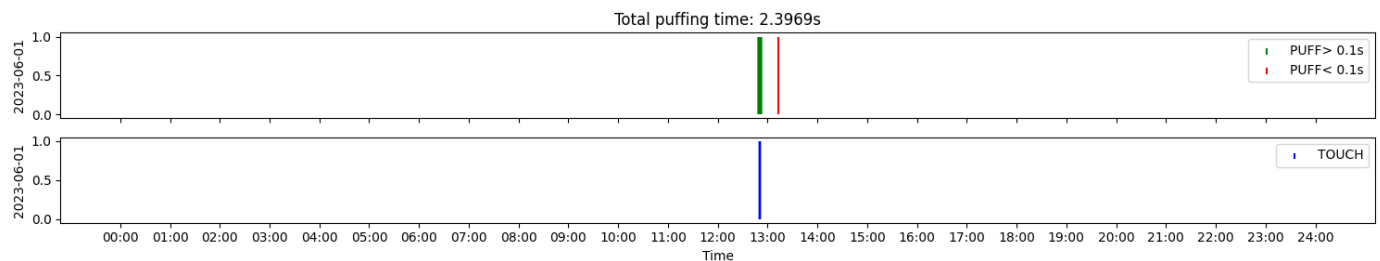

**Clear:** This button will clear all the text in RX monitor

**Save the plot:** For saving the plot, press the icon that is marked by red circle below.

Figure 1

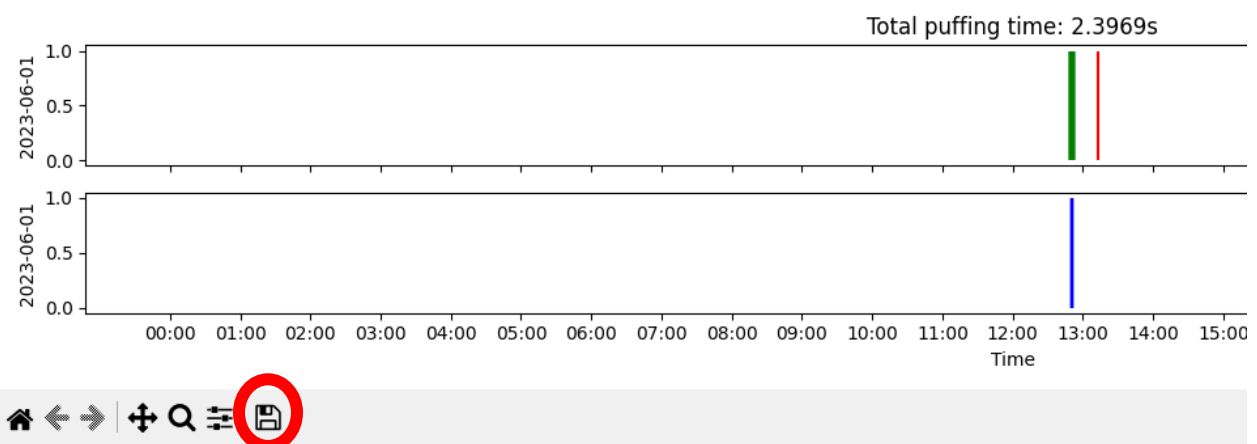

**N.B:** If any issue occurs while reading the data, please clear the Rx monitor by pressing clear button and try again. This may solve the issue.